

## OXYGEN NASAL CANNULAE/CATHETERS

**Medical device class; 2A**

### Intended use

Delivery of medical oxygen to patient

### Intended user


The device must only be used by clinical staff or healthcare professionals who have received training in oxygen therapy / respiratory procedures.

### Intended use environment

The devices are intended for use healthcare facilities and first response environments.

If the devices are to be used in homecare environments, they must be used under the supervision of a suitably trained healthcare professional.

The devices are single-use non-sterile

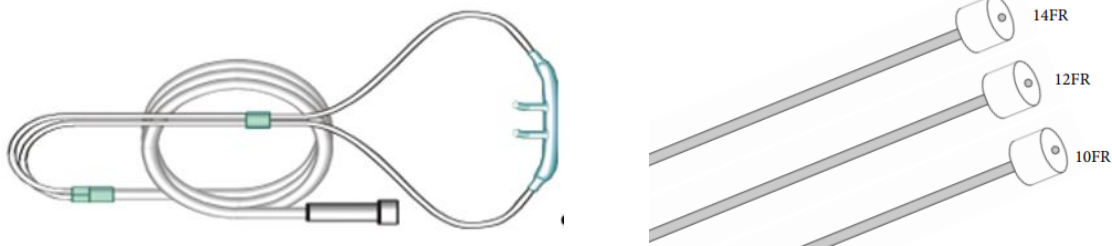

**GMDN; 35201**

### Material

- PVC DEHP FREE – NASAL CANNULAE
- PVC FREE – NASAL CATHETERS
- Where phthalates are in the materials the “phthalates” symbol is on the product label.
- Materials do not contain latex or substances derived from animals

### Instructions for use

NOTE. These directions are general guidelines intended for use by qualified medical personnel.

1. Attach oxygen supply tubing end to gas source and set prescribed flow.
2. Check for gas flow through the device
3. Place the nasal prongs in the patient's nose with the twin straps over the ears and around the neck.
4. Gently pull the slider on the strap up to the patient's chin until the device is secure.

### Storage

- ⇒ Do not store at extreme temperature or humidity.
- ⇒ Store in clean and dry place.
- ⇒ Do not use if the package open or damaged
- ⇒ Do not store in direct sunlight

### Cleaning

DISPOSABLE - FOR SINGLE PATIENT USE ONLY. NOT INTENDED FOR REPROCESSING

There is a risk of cross contamination if these devices are reused. They have not been validated for reprocessing or reuse.

### Warnings and cautions

1. Product may contain phthalates and could pose a risk for vulnerable patients including pregnant women or infants. Check the device label for the indication of the presence of phthalates.
2. Be sure all connections are secure.
3. Cleaning or reprocessing for reuse will raise the risk of cross contamination.
4. Meditech Systems Ltd will not take responsibility for use outside of the intended use or by other than the intended users.
5. FIRE RISK IF OXYGEN IS USED IN THE PRESENCE OF NAKED FLAME OR HEAT SOURCE LIKELY TO CAUSE IGNITION.

### Disposal

The devices may be disposed of using local procedures for contaminated medical waste.

The materials are capable of being recycled if a safe process for recycling medical waste is available.

### Ordering information:

Meditech Systems Ltd, Unit 3 Shrublands Estate, Shearstock, Shaftesbury, Dorset, SP7 9PT, England

[orders@meditechsystems.co.uk](mailto:orders@meditechsystems.co.uk) Tel: +44(0)1747 821546 Fax: +44(0)1747 825038

Facebook: [@meditechsystemslimited](https://www.facebook.com/meditechsystemslimited) LinkedIn: [meditech-systems-limited](https://www.linkedin.com/company/meditech-systems-limited) Twitter: [@MeditechLtd](https://twitter.com/MeditechLtd)

### EU Authorised Representative

















Clowrey Consultancy


Address; White House at Bridge, Holycross Village, Thurles, County Tipperary, E41 EC65, Eire.

Tel: +44(0)7888744362

E-mail: [clowreyconsult@gmail.com](mailto:clowreyconsult@gmail.com)

### Explanation of symbols:

	= Single Use. The definition of single use is for use on a single patient for a single procedure.		= Manufacturer
	= Do not reprocess / re-sterilise		= Temperature range
	= Use by date		= Consult instructions for use
	= Lot or batch number.		= Keep away from sunlight
	= Date of manufacture.		= Keep dry
	= Device catalogue number / product code / stock		= Contains latex
	= Caution		= Do not use of packaging is wet or damaged.
	= Device catalogue number / product code / stock		= Contains or Presence of Phthalates.


 = CE mark with notified body number

*This is an explanation of the symbols that may be found on Meditech Systems Ltd's products.*

### ADVERSE INCIDENT REPORTING

Adverse incidents may be reported to the Competent Authority in the country / region where the incident occurred.